

Introduction

Dear customer,

DIATOME Diamond Knives – Development, manufacturing and customer service since 1970

Our developments in this period:

ultra 45°	the diamond knife with a hydrophilic cutting edge, allowing high quality ultrathin sections of epoxi resin embedded samples			
ultra 35°	the low angle diamond knife for reduced compression and better structure preservation			
histo	the first diamond knife for cutting semithin sections for the observation in the optical microscope			
cryo	the diamond knife for cutting in the cryo-ultramicrotome			
Static Line Ionizer	for eliminating electrostatic charging			
cryo immuno	the knife with a large diamond platform fascilitating pick-up of cryo sections from sucrose infiltrated samples (Tokuyasu)			
CEMOVIS 25° CEMOVIS 35°	the low angle diamond knife for sectioning frozen hydrated samples (CEMOVIS)			
ultra AFM and cryo AFM	the knives for generating absolutely smooth and flat surfaces of biological and technical samples for the observation in the AFM			
ultra sonic	the oscillating diamond knife for eliminating compression and allowing best structure preservation			
accessories	such as pick-up loops, special forceps, special knife boats, etc.			

These developments were possible by the co-operation with you, our valued customers.

We are convinced that also in the future new developments are possible and that our knives may be adapted to the changing requirements.

Make use of our many years experience in perfecting our knives.

With a telephone call or an email we can inform you on any details you require.

We are looking forward to hear from you soon!

Your Diatome team

Content

Characteristics of Diatome Diamond Knives Resharpening and Exchange Service	4
Diatome Quality Guarantee Diatome Customer Service	5
DIATOME ultra knives	6 – 13
ultra 35°	
ultra semi	
ultra AFM	
ultra 35° Jumbo	
ultra sonic	
ultra 45°	
ultra 45° Jumbo	
trim 45 / trim 20 / trim 90	
DIATOME crue knives	14 10

DIATOME cryo knives

CEMOVIS 25° / CEMOVIS 35° cryo immuno cryo 35° (dry) cryo AFM cryo 35° (wet) cryo 45° (dry) cryo 45° (wet) trim 45 / trim 20 / trim 90

DIATOME histo knives

20-23

histo histo cryo (dry) histo cryo (wet) histo Jumbo

Characteristics of DiATOME Diamond Knives

Diatome knives are compatible with all ultramicrotomes

The boat is designed in such a manner that the water is horizontal when the clearance angle is set. This results in a stationary water surface and good reflection.

The shape of the boat offers you the advantage of an easy pick-up of the floating sections.

The surfaces of the Diatome diamond knives are hydrophilic. This property allows an easy wetting of the cutting edge, even with a very low water level (important for sectioning Lowicryls and other hydrophilic embedding media, as well as for water sensitive samples).

The diamond is high vacuum brazed onto a metal shank which is securely screwed into the boat.

The cementing material seals boat and knife. It is hydrophilic and aids the easy wetting of the cutting edge.

Resharpening and Exchange Service

Resharpened Diatome knives undergo the same stringent optical checking and sectioning test as new knives.

A diamond knife resharpened by Diatome is the same high quality as a new one!

We let you choose: You may have your knife resharpened, or exchanged against a new knife at a slightly elevated price.

In the exchange price a type change is possible (ultra to cryo, 45° to 35°, small knives to large ones, etc).

We guarantee a fast resharpening service for all our knives including the oldest knives ever purchased!

DiATOME Quality Guarantee

Before delivery, each knife, new or resharpened, is subject to extensive testing. We can therefore guarantee perfect sectioning (in the repective thickness range) over the entire cutting edge.

Should a knife not fulfil your high expectations, please kindly send it back with a short note of the problem.

DiATOME Customer Service

The great experience obtained during long years of development, manufacturing and use of diamond knives allows us to offer a unique sample testing service:

You send us biological or material research specimens of any nature.

We perform a sectioning test and send you the resulting ultrathin sections on grids, the surfaced sample, along with a report of how the results were obtained.

You may be assured that we treat your proprietary specimens with absolute discrection.

Please allow us to help you choose the appropriate knife type from our large range for your specific application.

DiATOME ultra knives and their applications

- Highest quality diamonds and optimal crystal orientation guarantee perfect ultrathin sections and a durable cutting edge
- Section pick-up is facilitated as the boat is horizontal allowing the water to completely fill the boat all the way round
- A hydrophilic surface makes it easy to wet the cutting edge, even with low water level

Knife type	Knife angle	Size [mm]	Thickness range [nm]	Boat type	Code	Application
ultra 35°	35°	1.5 2.0 2.5 3.0 3.5 4.0	30-200	Standard boat	DU3515 DU3520 DU3525 DU3530 DU3535 DU3540	Biological and materials science specimens
ultra semi	35°	3.0	50 – 500	Standard boat	DU3530-semi	 Alternating sectioning ultrathin semithin
ultra AFM	35°	2.0 3.0	15 – 100	Small boat	DU3520-AFM DU3530-AFM	Sample surfacing for AFM
ultra 35° Jumbo	35°	3.0	50 – 200	Jumbo boat	DUJ3520 DUJ3530	 Section series for 3D reconstruction, STEM
ultra sonic	35°	3.0	15 – 100	Special boat	DUS3530	 Compression free sections Best structure preservation Biological and materials science specimens
ultra 45°	45°	1.5 2.0 2.5 3.0 3.5 4.0	30-200	Standard boat	DU4515 DU4520 DU4525 DU4530 DU4535 DU4540	Routine sectioning of biological and materials science specimens
ultra 45° Jumbo	45°	3.0	50-200	Jumbo boat	DUJ4530	Section series for 3D reconstruction, STEM
trim 45 trim 20 trim 90	45°			Triangular holder	DTB45 DTB20 DTB90	Trimming biological and materials science specimens

ultra

ultra 35°

In 1989 considerably reduced compression, smoother section surfaces and improved structural preservation thanks to the use of our ultra 35° knives was demonstrated (J. C. Jésior, Scanning Microscopy Supplement 3, pp. 17 – 153, 1998).


In the meantime, a large number of researchers have recognized the advantages of 35° knives, in particular for sectioning biological specimens of all kind, non-homogenous specimens, non decalcified bone, dental material, etc.

The ultra 35° knives are perfect for sectioning relatively soft materials research specimens including metals and polymers, as well as mixed specimens such as polymers filled with nanoparticles, brittle materials such as catalysts, crystals, semiconductors, etc (G. Mahon et al., Microscopy Research and Technique, Vol. 31, pp. 267 – 274, 1995, S. R. Glanvill, Microscopy Research and Technique, Vol. 31, pp. 275 – 284, 1995, P. Swab et al., Mat. Res. Soc. Symp. Proc. Vol. 115, pp. 229 – 234, 1989, P. Schubert-Bischoff et al., Microscopy and Microanalysis, proceedings, page 359, 1997).


The ultra 35° knife has demonstrated it's usfulness as a standard knife for the majority of applications in both biological and materials research.



EM micrograph of an ultramicrotomed section of the anodic alumina film formed on Al-2 wt%Cu alloy. Scale bar = 100 nm. Xiarong Zhou, School of Materials, University of Manchester.

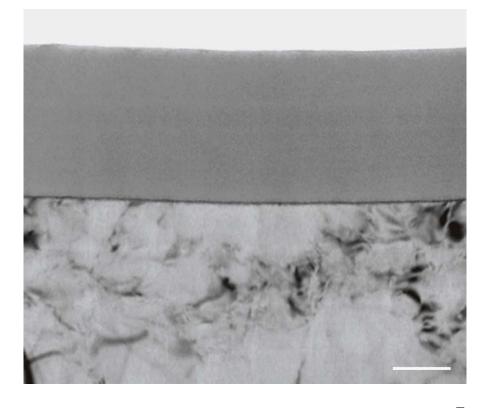

▼
Rat brain x 18'000
Werner Graber, Anatomisches Institut, Bern.

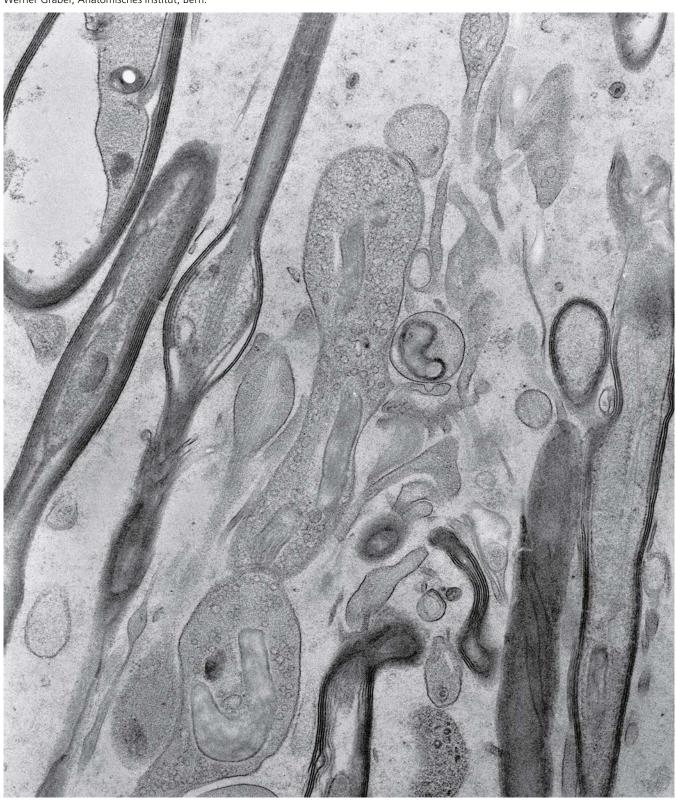


ultra AFM

In order to achieve the best results for AFM investigation, only the highest quality diamond knives should be used (P. H. Vallotton et al., J. Biomater. Sci. Polymer Edn., Vol. 6, No. 7, pp. 609 – 620, 1994. N. Matzko et al., Journal of Structural Biology 146, pp. 334 – 343, 2004).

Our ultra AFM knives are made of highest quality to ensure the increased quality requirements of AFM investigation. They preduce extremely smooth sample surfaces and guarantee the best possible structure preservation.



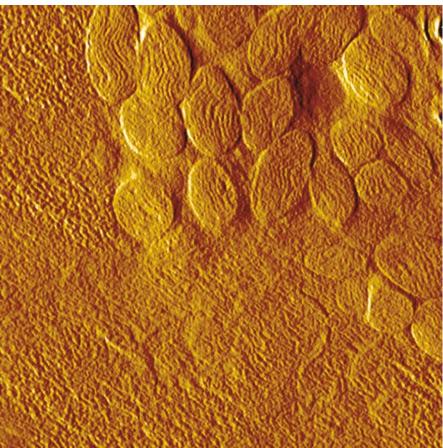

AFM amplitude image of the muscle of cat's mite Otodectes cynotis. The contrast covers amplitude variation in the 1-3 nm range. Size of the whole image equals 4.6 microns. Nadejda Borisovna Matsko, Institut für angewandte Physik, ETH Zürich.

ultra sonic

- Thinner sections
- No compression
- Best structure preservation

The patented ultra sonic knife allows the cutting of ultrathin sections free of compression (D. Studer et al., Journal of Microscopy, Vol. 197, pp. 94 – 100, 2000).


With correct setting of frequency and amplitude the sections become as long as the height of the sample.

Best results are not only achieved with biological samples, but also with polymers (J. S. Vastenhout et al., Microscopy and Microanalysis 8, 2001. J. S. Vastenhout et al., Microscopy Today, pp. 20–21, 2006).

We have tested the ultra sonic knife with the following samples:

- Biological samples in Epon, Araldite, EM Bed, etc.
- Biological samples in acrylic resins (Lowicryl, LR White).

 Rigid polymers such as PS, PMMA, ABS, HIPS, modified PP, blends of various kind.

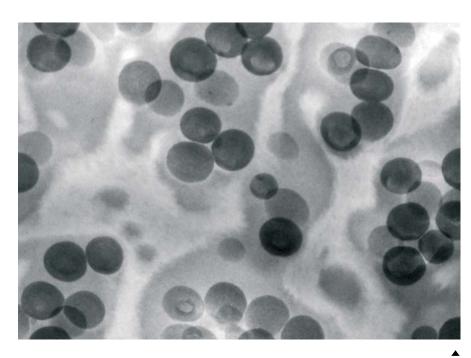



Neuropil of a HP frozen / freeze substituted hippocampus slice culture, Daniel Studer and Werner Graber, Inst. of Anatomy, University of Bern. Sectioned with the ultra sonic knife.




ultra 45°

Acknowledged as the appropriate knife angle for routine sectioning of both biological and materials research specimens, it represents a balanced compromise between section quality and durability.

For the sectioning of a number of hard materials such as ceramics, semiconductors, oxides etc, with the use of the ultra 45° knife a longer service time may be expected. Kindly contact us and make use of our long years experience in all ultramicrotomy applications.



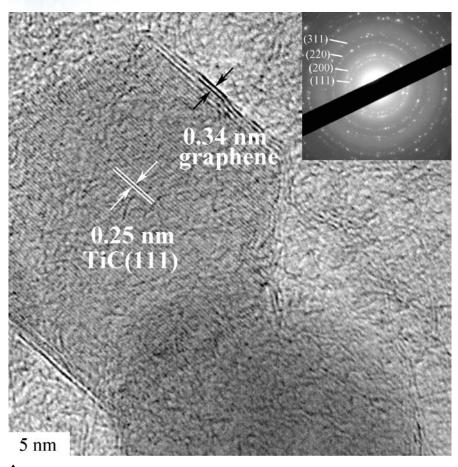

HRTEM image of a TiC microsphere showing nanocrystalline carbon coating on the TiC grain (confirmed by the electron diffration pattern).
J. Zhong et al, Journal of the European Ceramic Society 32, pp. 3407 – 3414, 2012.

trim 45, trim 20 and trim 90

- Rapid and precise trimming
- Shiny block faces and pyramidal sides
- Sample surface aligned with the cutting direction
- Durable cutting edge

For successful ultramicrotomy in biology and materials science, precise trimming is mandatory.


The Diatome knives trim 45, trim 20 and trim 90 will fulfil your trimming requirements, allowing quick, easy and accurate trimming.

A well trimmed sample is a precondition for perfect section ribbons.

Trimming with the trim 45, trim 20 and trim 90 improves the service time of your diamond knives.

When a sample is trimmed with glass knives or razor blades hard particles may be introduced in the surface or the pyramidal sides. These particles immediately cause knife marks in the diamond knives.



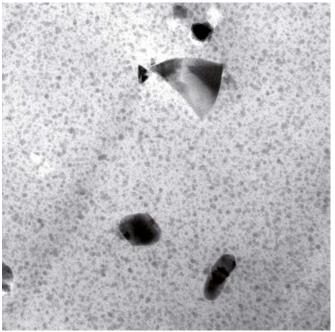

SiC, SiO2, TiO2 and AlO3 nanoparticles in polymer matrix. Claudia Mayrhofer, TU Graz.

→ 0.5 μm

DiATOME cryo knives and their applications

- Thinner cryo sections
- Perfect cryosections from ultrathin to semi with the same knife
- Minimal compression and best structure preservation
- Highest quality diamonds and optimal crystal orientation guarantee perfect ultrathin sections and a durable cutting edge


Knife type	Knife angle	Size [mm]	Thickness range [nm]	Boat type	Code	Application
CEMOVIS 25°	25°	3.0	30 – 150	Triangular holder	DCEM2530	 Frozen hydrated samples (CEMOVIS)
CEMOVIS 35°	35°	3.0	30 – 150	Triangular holder	DCEM3530	 Frozen hydrated samples (CEMOVIS)
cryo immuno	35°	2.0 3.0	30-300	Support triangulaire	DCIMM3520 DCIMM3530	 Sucrose infiltrated samples (Tokuyasu)
cryo 35° (dry)	35°	1.5 2.0 2.5 3.0 3.5 4.0	30 – 300	Triangular holder	DCO3515 DCO3520 DCO3525 DCO3530 DCO3535 DCO3540	• Polymers, rubber, paints, etc.
cryo AFM	35°	2.0 3.0	20 – 100	Triangular holder	DCO3520-AFM DCO3530-AFM	 Sample planing for AFM imaging
cryo 35° (wet)	35°	1.5 2.0 2.5 3.0 3.5 4.0	30 – 300	Small cryo boat	DCM3515 DCM3520 DCM3525 DCM3530 DCM3535 DCM3540	 Wet sectioning of polymers with DMSO/ water mixture
cryo 45° (dry)	45°	1.5 2.0 2.5 3.0 3.5 4.0	30 – 300	Triangular holder	DCO4515 DCO4520 DCO4525 DCO4530 DCO4535 DCO4540	 Routine dry cryo sectioning of polymers
cryo 45° (wet)	45°	1.5 2.0 2.5 3.0 3.5 4.0	30 – 300	Small cryo boat	DCM4515 DCM4520 DCM4525 DCM4530 DCM4535 DCM4540	 Routine wet cryo sectioning of polymers with DMSO/water
trim 45 trim 20 trim 90	45°			Triangular holder	DTB45 DTB20 DTB90	 Trimming biological and materials science samples

CEMOVIS 25° CEMOVIS 35°

The CEMOVIS 35° knife and the CEMOVIS 25° knife are designed for sectioning frozen hydrated specimens. The 25° angle results in the least possible compression and the best structure preservation (H.M. Han et al., Journal of Microscopy, Vol. 230, Pt. 2, pp. 167–171, 2007).

Please note: best results are achieved at low humidiy, when the cryo-ultramicrotome is placed in a glovebox and the sections attached by electrostatic force (J. Pierson et al., Journal of Structural Biology 169, pp. 219 – 225, 2010).



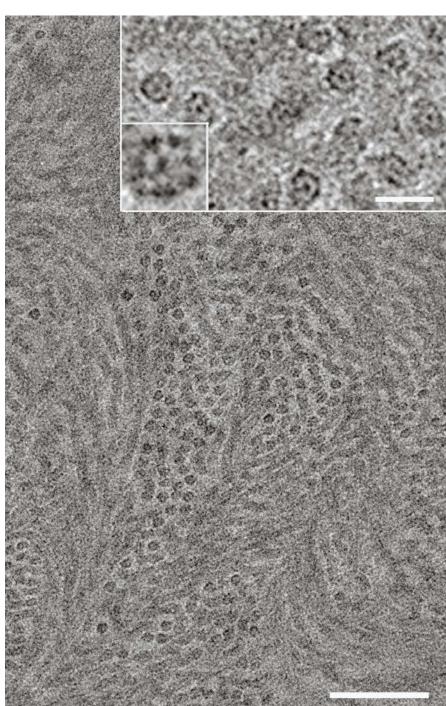

▲ High resolution electron micrograph of vitreous section of keratin intermediate filaments in the midportion of stratum corneum of human epidermis.

The fine structure of the keratin filaments is well resolved and their molecular organisation is seen in favourable cases (inset).

Scale bar = 100 nm.

Scale bar = 100 nm.
Scale bar inset = 20 nm.
Ashraf Al-Amoudi, Laboratoire d'Analyse
Ultrastructurale, Lausanne.

cryo immuno


The first cryo knife with a diamond platform, guarantees the best possible sectioning for sucrose infiltrated samples (Tokuyasu).

The diamond platform guarantees an easy and gentle section pick-up.


The sections are collected directly from the diamond surface using a loop and a sucrose/methyl-cellulose droplet (W. Liou et al., Histochemistry and Cell Biology, Vol. 106, pp. 41 – 55, 1996. P. J. Peters et al., Current Protocols in Cell Biology, pp. 4.7.1 – 4.7.18, 2006).

The 35° angle leads to a considerable reduction in mechanical stresses and therefore to improved structure preservation in sucrose infiltrated samples (E. Bos et al., Journal of Structural Biology 175, pp. 62 – 72, 2011).



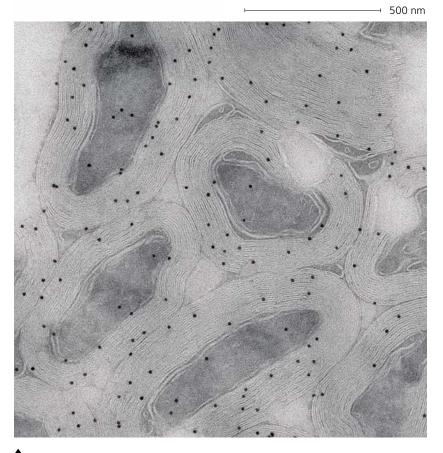

Mouse optic nerve, immunolabeling of the major myelin protein proteolipid protein (PLP), 10 nm gold. Wiebke Möbius, Dept. of Neurogenetics, EM Core Facility, MPI of Experimental Medicine, Göttingen.

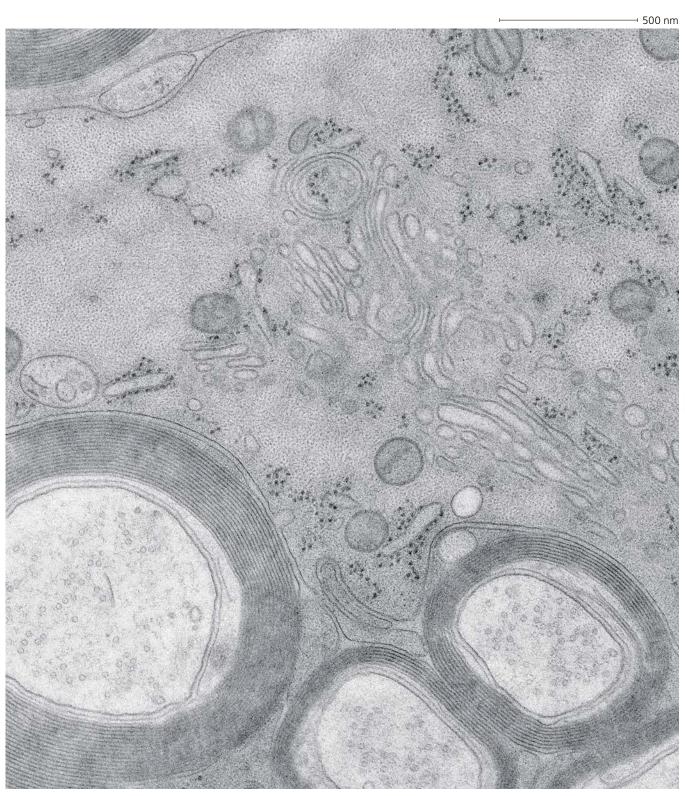

Ultrastructure of an astrocyte in the mouse optic nerve, high-pressure frozen. In the astrocyte the glial fibrillary acidic protein is visible.
Wiebke Moebius, Dept. of Neurogenetics, EM Core Facility, MPI of Experimental Medicine, Göttingen.

cryo 35° cryo 45°

The cryo 35° knife has demonstrated its usefulness as a standard knife for the low temperature sectioning of polymers, rubber, paints, etc.

The 35° angle represents a balanced compromise between section quality and durability.

The cryo 35° and cryo 45° knife mounted in the triangular holder is suitable for dry cryosectioning.

The cryo 35° and cryo 45° knife mounted in the trough are used for sectioning with fluids such as a DMSO/water mixture.

The cryo 45° knife is well suited for routine cryo sectioning.




▲ Styrene-butadiene block copolymer x 25′000 Ronald Walter, BASF Aktiengesellschaft, Polymer Physics, D-67056 Ludwigshafen.



cryo AFM

Our cryo AFM knives are made of hightest quality to ensure the increased quality requirements of AFM investigation. They preduce extremely smooth sample surfaces and guarantee the best possible structure preservation.



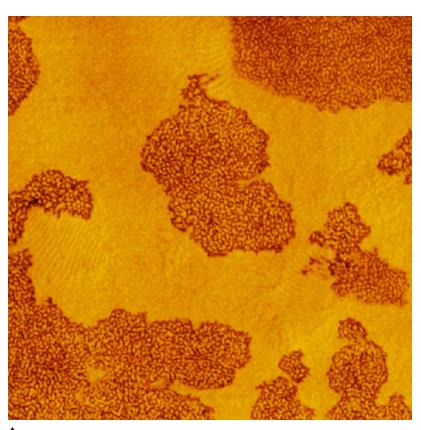

Morphology of a blend of two SBS block copolymers with different chain-architecture. AFM tapping mode, phase image, image size = 3 x 3 µm. Rameshwar Adhikari, Institut für Werkstoffwissenschaft, Martin-Luther-Universität, Halle-Wittenberg.

trim 45, trim 20 and trim 90

For successful ultramicrotomy in biology and materials science, precise trimming is mandatory.

The Diatome knives trim 45, trim 20 and trim 90 diamond blades will fulfil all your trimming requirements, allowing quick, easy and accurate trimming.

Due to the extreme sharpness of our diamond blades, less mechancial damage is applied to the sample during trimming. Very shiny sample faces and precise sides are the result.



DIATOME histo knives and their applications

- High quality diamonds guarantee perfect sections and a durable cutting edge
- Easy wetting cutting edge
- Optimised serial sectioning
- Large boat for easy pick-up


Knife type	Knife angle	Size [mm]	Thickness range [µm]	Boat type	Code	Application
histo	45°	4.0 6.0 8.0	0.2-2	Standard histo boat	DH4540 DH4560 DH4580	 Sectioning biological and industrial materials specimens for optical microscopy
histo Jumbo	45°	6.0 8.0	0.2-2	histo Jumbo boat	DHJ4560 DHJ4580	• Large boat for serial sectioning
histo cryo (dry)	45°	4.0 6.0 8.0	0.2-2	Triangular histo holder	DHCO4540 DHCO4560 DHCO4580	Sectioning biological and materials specimens for optical microscopy
histo cryo (wet)	45°	4.0 6.0 8.0	0.2-2	Small histo boat	DHCM4540 DHCM4560 DHCM4580	Wet sectioning of polymers (with DMSO/water mixture)

histo

The histo knife is designed for the sectioning of hard and soft biological and material research specimens, non embedded or embedded in acrylic or epoxy resins (O. L. Reymond, Bas. Appl. Histochem. 30, pp. 487 – 494, 1986).

Our histo knives can be used on all ultramicrotomes as well as on microtomes with a retraction of the specimen arm in the return phase.




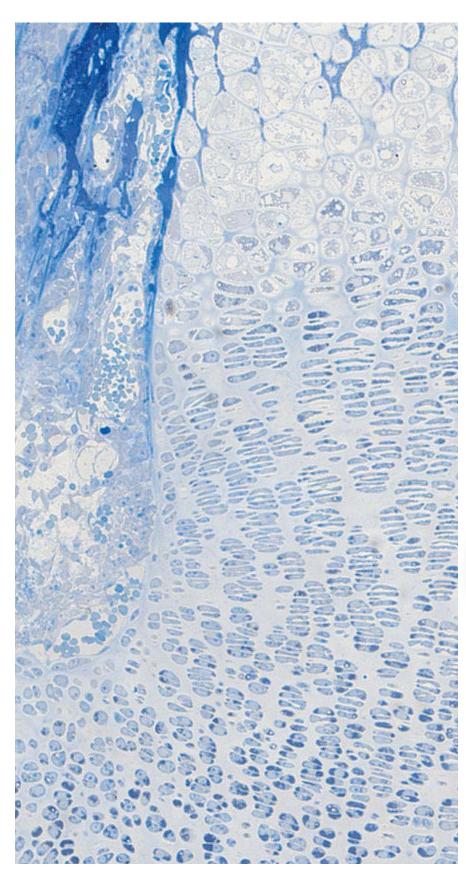

► Nondecalcified rat bone. Scale: 35 mm = 100 μm. Daniel Studer, Anatomisches Institut, Bern.

histo cryo

Our histo cryo knives are delivered with a boat for wet sectioning using a DMSO/ water mixture or in triangular holders for dry sectioning.

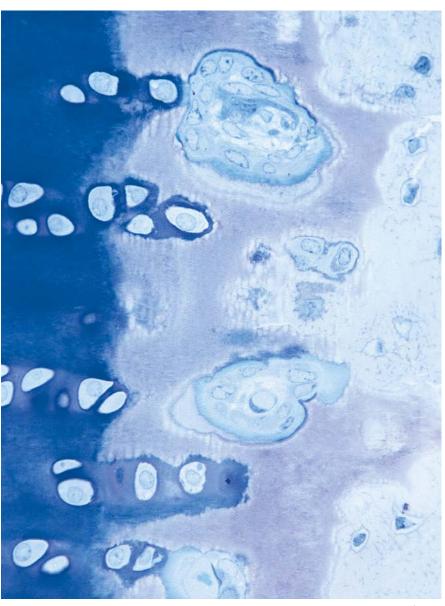





▲ E. B. Hunziker, M. E. Müller, Institute for Biomechanics, University of Bern. Rabbit joint, calcified cartilage/bone. X760

histo Jumbo


For 3D reconstruction it is imperative not to loose a single section (M. J. F. Blumer et al., Journal of Neuroscience Methods 120, pp. 11 – 16, 2002).

The large Jumbo boat as well as the adhesive (Pattex compact by Henkel) applied to the upper side of the sample block increase the distinct advantages of our histo knives.

They allow:

- easy production of section ribbons (0.5 – 2 µm)
- no section loss
- no folding
- the same orientation of all sections
- easy collection of section ribbons
- multiple ribbons on one glass slide
- perfect for immuno-histo-chemistry.



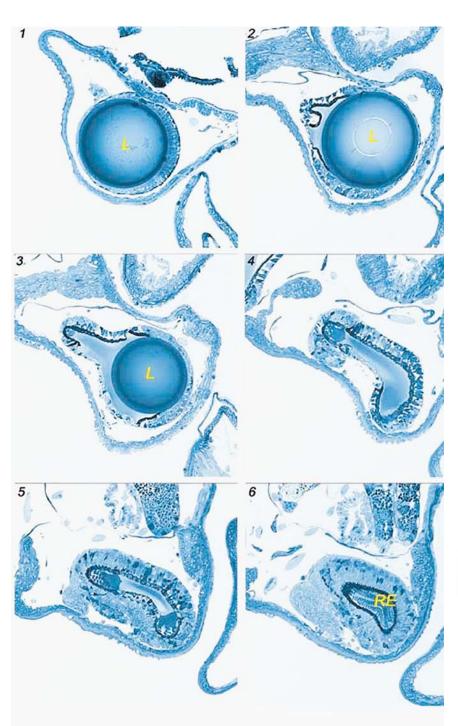

 \blacktriangle Eye of A. peroni: photographic sequence of some semithin sections of a complete serie through the eye. L = lens; RE = retina.

Michael J. F. Blumer, Institut für Zoologie, Universität Wien.

Reprinted from: Ribbons of semithin sections: an advanced method with a new type of diamond knife. Journal of Neuroscience Methods 120, pp. 11 – 16, 2002, with permission from Elsevier.